




# iCell<sup>®</sup> Cardiomyocytes<sup>2</sup> User's Guide



## Trademarks

iCell is a registered trademark, and Cellular Dynamics and the  logo are trademarks of Cellular Dynamics International, Inc. and may not be used without the express written permission of Cellular Dynamics International, Inc. (CDI).

All other brands, product names, company names, trademarks, and service marks are the properties of their respective owners.

## Restrictions and Liabilities

This document is provided "as is." CDI assumes no responsibility for any typographical, technical, or other inaccuracies in this document. CDI reserves the right to periodically change information that is contained in this document; however, CDI makes no commitment to provide any such changes, updates, enhancements, or other additions to this document to you in a timely manner or at all.

THE ICELL CARDIOMYOCYTES<sup>2</sup> ARE PROVIDED "AS IS" TO YOU. TO THE FULLEST EXTENT PERMITTED BY APPLICABLE LAW, CDI MAKES NO, AND DISCLAIMS ALL, REPRESENTATIONS, WARRANTIES, CONDITIONS OR COVENANTS, EXPRESS OR IMPLIED, WITH RESPECT TO ANY PRODUCT REFERENCED HEREIN OR PERFORMANCE OF ANY SERVICES REFERENCED HEREIN, INCLUDING BUT NOT LIMITED TO, ANY EXPRESS OR IMPLIED WARRANTIES OR CONDITIONS OF FITNESS FOR A PARTICULAR PURPOSE, NON-INFRINGEMENT, MERCHANTABILITY, DURABILITY, TITLE, OR RELATED TO THE PERFORMANCE OR NON-PERFORMANCE OF ANY PRODUCT REFERENCED HEREIN OR PERFORMANCE OF ANY SERVICES REFERENCED HEREIN).

This document might contain references to third-party sources of information, hardware or software, products or services and/or third-party websites (collectively the "Third-Party Information"). CDI does not control, and is not responsible for, any Third-Party Information, including, without limitation the content, accuracy, copyright compliance, compatibility, performance, trustworthiness, legality, decency, links, or any other aspect of Third-Party Information. The inclusion of Third-Party Information in this document does not imply endorsement by CDI of the Third-Party Information or the third party in any way.

CDI does not in any way guarantee or represent that you will obtain satisfactory results from using iCell Cardiomyocytes<sup>2</sup> as described herein. You assume all risk in connection with your use of iCell Cardiomyocytes<sup>2</sup>.

## Conditions of Use

iCell Cardiomyocytes<sup>2</sup> are for life science research use only and subject to the use restrictions as contained in Appendix A. You are responsible for understanding and performing the protocols described within. CDI does not guarantee any results you may achieve. These protocols are provided as CDI's recommendations based on its use and experience with iCell Cardiomyocytes<sup>2</sup>.

## Origin

iCell Cardiomyocytes<sup>2</sup> are manufactured in the United States of America.

## Copyright Notice

© 2016 Cellular Dynamics International, Inc. All rights reserved. This document may not be reproduced, distributed, modified or publicly displayed without the express written permission of Cellular Dynamics International, Inc.

## Revision History

Version 0.4: August 2016

Version 0.3: September 2015

Version 0.2: July 2015

Version 0.1: March 2015

## Table of Contents

<b>Before You Begin</b> .....	<b>ii</b>
<b>Chapter 1. Introduction</b> .....	<b>1</b>
Components Supplied by Cellular Dynamics .....	2
Required Equipment and Consumables.....	2
Technical Support and Training.....	3
Workflow Diagram .....	4
<b>Chapter 2. Handling and Storage</b> .....	<b>5</b>
Handling iCell Cardiomyocytes <sup>2</sup> .....	5
Handling iCell Cardiomyocytes Media.....	5
<b>Chapter 3. Preparing Cell Culture Surfaces</b> .....	<b>6</b>
Preparing 0.1% (w/v) Gelatin Solution.....	6
Preparing the Gelatin Cell Culture Vessel.....	6
<b>Chapter 4. Thawing Media</b> .....	<b>7</b>
<b>Chapter 5. Thawing iCell Cardiomyocytes<sup>2</sup></b> .....	<b>8</b>
<b>Chapter 6. Plating iCell Cardiomyocytes<sup>2</sup></b> .....	<b>10</b>
<b>Chapter 7. Maintaining iCell Cardiomyocytes<sup>2</sup></b> .....	<b>11</b>
<b>Appendices</b> .....	<b>12</b>
Appendix A. Intellectual Property Rights, Use Restrictions, and Limited License...12	
Appendix B. Product Provided “AS IS” .....	12
Appendix C. Limited Liability.....	13

## Before You Begin

Notes


- Immediately transfer the frozen vials to liquid nitrogen storage.
- Read this entire iCell® Cardiomyocytes<sup>2</sup> User's Guide before handling or using iCell Cardiomyocytes<sup>2</sup>.
- iCell Cardiomyocytes<sup>2</sup> are for life science research use only. See Appendix A for more information and other restrictions.
- A Safety Data Sheet (SDS) for dimethyl sulfoxide (DMSO), in which iCell Cardiomyocytes<sup>2</sup> are frozen, is available online at [www.cellulardynamics.com/lit/](http://www.cellulardynamics.com/lit/) or on request from Cellular Dynamics International. Only technically qualified individuals experienced in handling DMSO and human biological materials should access, use, or handle iCell Cardiomyocytes<sup>2</sup>.

## Chapter 1. Introduction

Cellular Dynamics International's (CDI) iCell Cardiomyocytes<sup>2</sup> are highly purified, human cardiomyocytes derived from induced pluripotent stem (iPS) cells using CDI's proprietary differentiation and purification protocols. iCell Cardiomyocytes<sup>2</sup> are a mixture of spontaneously electrically active atrial-, nodal-, and ventricular-like myocytes with typical biochemical, electrophysiological, and mechanical characteristics and expected responses upon exposure to exogenous agents. Thus, these cells provide a reliable source of human cardiomyocytes suitable for use in targeted drug discovery, toxicity testing, and other life science research.

When thawed and plated with iCell Cardiomyocytes Plating Medium and maintained in iCell Cardiomyocytes Maintenance Medium as instructed in this User's Guide, iCell Cardiomyocytes<sup>2</sup> will begin to beat spontaneously within 24 hours. When seeded at appropriate densities, iCell Cardiomyocytes<sup>2</sup> also will form electrically connected syncytial layers that beat in synchrony and can be used for electrophysiology interrogation via such techniques as microelectrode array and impedance measurements on day 4 post- thawing.

iCell Cardiomyocytes Maintenance Medium is antibiotic-free and has been specially formulated to maintain the health and function of the cardiomyocytes while limiting the proliferation of the small percentage of non-cardiomyocyte cells. iCell Cardiomyocytes<sup>2</sup> therefore can be maintained in culture for at least 14 days in the Maintenance Medium without appreciable loss of purity, enabling longer term studies. Thus, the combination of CDI's purification process and adherence to the procedures described in this User's Guide makes additional use of antibiotics unnecessary.

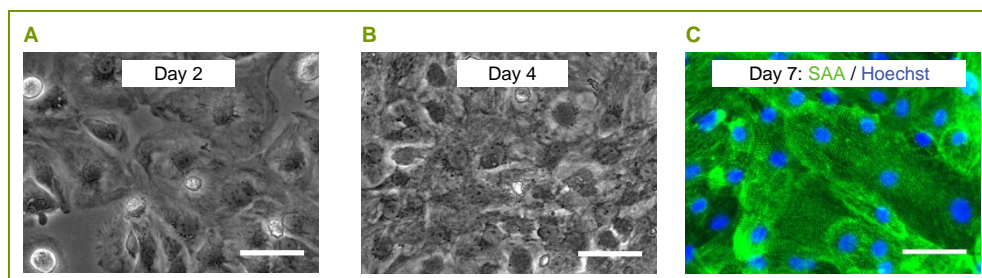

**Figure 1: iCell Cardiomyocytes<sup>2</sup> Represent a Highly Pure Population of Human Myocytes**

*These images show iCell Cardiomyocytes<sup>2</sup> at days 2, 4 and 7 post-plating. iCell Cardiomyocytes<sup>2</sup> form a syncytial layer in culture within 4 days post-plating and show the appropriate cardiac markers as demonstrated by immunocytochemistry: sarcomeric alpha-actinin (SAA, green) and Hoechst (blue). Scale bar = 50 µm.*

## Components Supplied by Cellular Dynamics

Item	Catalog Number
iCell Cardiomyocytes <sup>2†</sup>	CMC-100-012-000.5 (0.5 unit) CMC-100-012-001 (1 unit)
iCell Cardiomyocytes Plating Medium <sup>†</sup>	CMM-100-110-001 (1 unit) CMM-100-110-005 (5 unit)
iCell Cardiomyocytes Maintenance Medium <sup>†</sup>	CMM-100-120-001 (1 unit) CMM-100-120-005 (5 unit)
iCell Cardiomyocytes <sup>2</sup> User's Guide <sup>†</sup>	
Certificate of Testing <sup>#</sup>	
Certificate of Origin If required for shipping purposes	

<sup>†</sup> Safety Data Sheets and User's Guide available online at [www.cellulardynamics.com/lit/](http://www.cellulardynamics.com/lit/)

<sup>#</sup> Available online at [www.cellulardynamics.com/cot/](http://www.cellulardynamics.com/cot/)

## Required Equipment and Consumables

Item	Vendor	Catalog Number
<b>Equipment</b>		
37°C Water Bath	Multiple Vendors	
Biological Safety Cabinet with UV Lamp	Multiple Vendors	
Cell Culture Incubator	Multiple Vendors	
Hemocytometer or Automated Cell Counter*	Multiple Vendors	
Liquid Nitrogen Storage Unit	Multiple Vendors	
Phase Contrast Microscope	Multiple Vendors	
Pipettors	Multiple Vendors	
Tabletop Centrifuge	Multiple Vendors	
<b>Consumables</b>		
50 ml Centrifuge Tubes	Multiple Vendors	
Cell Culture Vessels, Sterile, TC Grade	Multiple Vendors	
Gelatin, Type A, Cell Culture Tested	Multiple Vendors	
Pipettes	Multiple Vendors	
Sterile Distilled Water	Multiple Vendors	
Trypan Blue	Gibco	15250

\* Ensure automated cell counter is appropriately calibrated before use.

---

## Technical Support and Training

CDI's Technical Support Scientists have the necessary laboratory and analytical experience to respond to your inquiries. In-lab training may be available upon request.


**Telephone** (877) 320-6688 (US toll-free) / (608) 310-5100 x5  
Monday - Friday, 8:30 am - 5:00 pm US Central Time

**Fax** (608) 310-5101

**Email** [support@cellulardynamics.com](mailto:support@cellulardynamics.com)

## Workflow Diagram

Notes

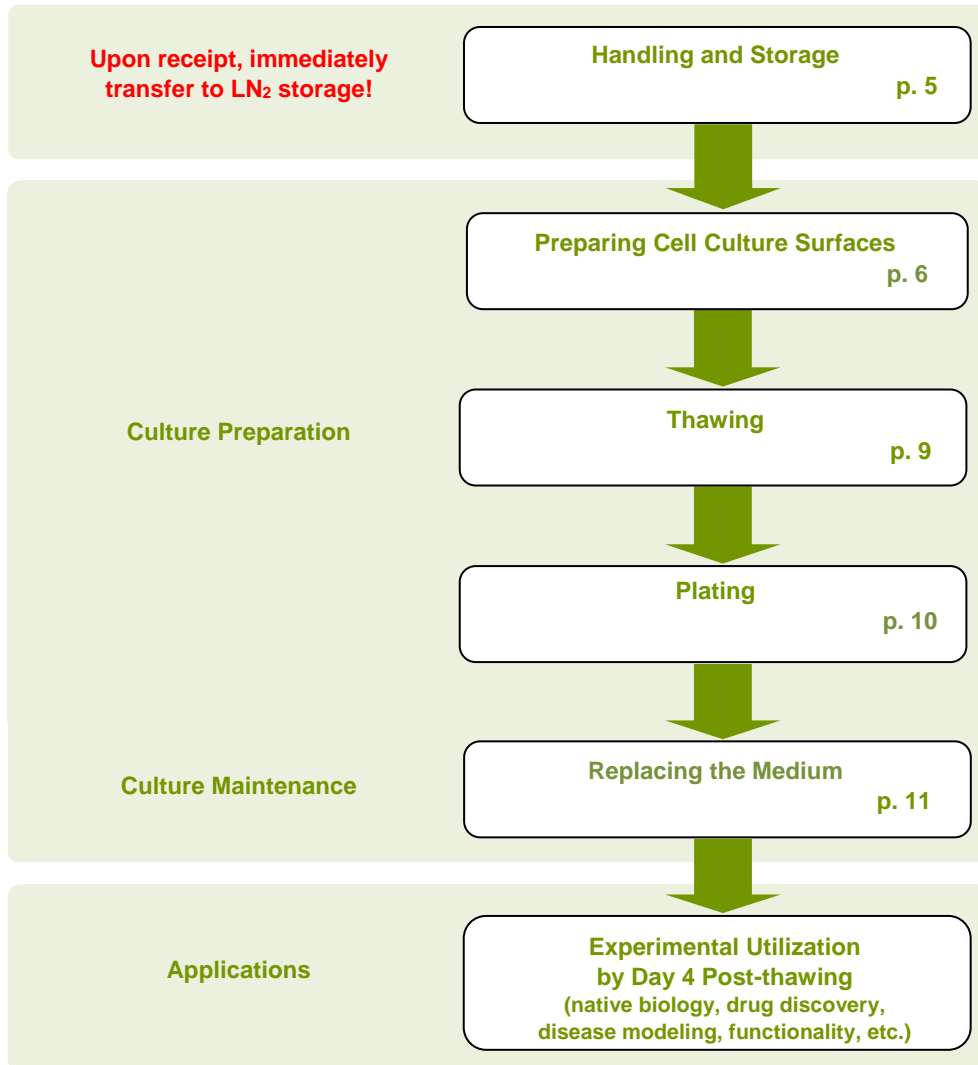



## Chapter 2. Handling and Storage

---

### Handling iCell Cardiomyocytes<sup>2</sup>

iCell Cardiomyocytes<sup>2</sup> are provided as cryopreserved single-cell suspensions in 1.5 ml cryovials. Upon receipt, directly transfer the cryobox containing iCell Cardiomyocytes<sup>2</sup> to the vapor phase of a liquid nitrogen storage dewar. CDI strongly recommends transferring the entire cryobox into the storage rack to avoid transferring individual vials.



*It is critical to maintain cryopreserved iCell Cardiomyocytes<sup>2</sup> at a stable temperature. Minimize exposure of cryopreserved iCell Cardiomyocytes<sup>2</sup> to ambient temperature when transferring vials to liquid nitrogen storage.*

---

### Handling iCell Cardiomyocytes Media

iCell Cardiomyocytes Plating Medium and iCell Cardiomyocytes Maintenance Medium are shipped frozen on dry ice. Upon receipt, store iCell Cardiomyocytes media at -20°C until ready for use and at 4°C for up to 2 weeks post-thawing. If media will be used for longer than 2 weeks, aliquot and freeze again after the initial thaw. Do not subject media to more than a single refreeze and thaw cycle.

## Chapter 3. Preparing Cell Culture Surfaces

iCell Cardiomyocytes<sup>2</sup> will plate and function on a variety of substrates including freshly prepared gelatin and fibronectin, which have been shown to support attachment, viability, and function of iCell Cardiomyocytes<sup>2</sup> with similar efficiencies.

CDI provides application protocols and notes that recommend assay-specific substrates. See [www.cellulardynamics.com/lit/](http://www.cellulardynamics.com/lit/) for a list of available application protocols for iCell Cardiomyocytes<sup>2</sup>. Regardless of the substrate of choice, prepare plating surfaces before thawing iCell Cardiomyocytes<sup>2</sup>.

### Preparing 0.1% (w/v) Gelatin Solution

1. Prepare bottles for mixing the 0.1% gelatin solution.
  - Do not use bottles that have previously contained or been exposed to detergent.
  - Acid-wash new bottles before use.
  - Use only high-purity, 18 M $\Omega$  sterile water for washing bottles.
2. Prepare 0.1% gelatin solution by dissolving 500 mg of gelatin (Type A, powder, cell culture tested) in 500 ml of sterile water in each bottle.
3. Sterilize the gelatin solution by autoclaving using liquid cycle. The sterile 0.1% gelatin solution is stable for 3 months when prepared as described above and stored at room temperature.

**Note:** Alternatively, use commercially available 0.1% gelatin in water (Stem Cell Technologies, Cat. No. 07903).

### Preparing the Gelatin Cell Culture Vessel

1. Select the cell culture vessel appropriate for your experimental use. Add the volume of 0.1% gelatin solution specified in the table below. Scale volumes appropriately for other vessel formats.

Culture Vessel	Surface Area (cm <sup>2</sup> )	Volume of 0.1% Gelatin Solution (ml)
12-well Cell Culture Plate	3.8	1.2
96-well Cell Culture Plate	0.32	0.1

**Table 1: Summary of Useful Volumes and Measures**

All volumes and measures are **per well**, if applicable.

**Note:** For glass coverslips for immunocytochemistry or electrophysiological applications, see the iCell Cardiomyocytes Application Protocols available online at [www.cellulardynamics.com/lit/](http://www.cellulardynamics.com/lit/).

2. Incubate the vessel(s) in a 37°C cell culture incubator for at least 1 hour.
3. Aspirate gelatin solution immediately before addition of the cell suspension.



Do not allow the gelatin-coated surface to dry.

## Chapter 4. Thawing Media

iCell Cardiomyocytes Plating Medium (Plating Medium) and iCell Cardiomyocytes Maintenance Medium (Maintenance Medium) have been specially formulated to maximize the cell viability and recovery at thaw, and to maintain the health and function of iCell Cardiomyocytes<sup>2</sup> in culture over time, respectively. Thaw and store the media as follows:

1. 24 hours before use, thaw the media overnight at 4°C.
2. Prepare aliquots of media and store at 4°C for up to 2 weeks.

**Note:** *The medium aliquots can be stored at -20°C. Do not thaw and refreeze the medium aliquots multiple times.*

## Chapter 5. Thawing iCell Cardiomyocytes<sup>2</sup>

Maintain iCell Cardiomyocytes<sup>2</sup> in liquid nitrogen until immediately before thawing to ensure maximal performance of the cells. Complete the following steps of the thawing procedure in a time-efficient manner to facilitate optimal iCell Cardiomyocytes<sup>2</sup> viability and performance.

**Note:** Thaw no more than 3 vials of iCell Cardiomyocytes<sup>2</sup> at one time.

1. Equilibrate the Plating Medium at room temperature before thawing iCell Cardiomyocytes<sup>2</sup>.
2. Remove the iCell Cardiomyocytes<sup>2</sup> cryovial from the liquid nitrogen storage tank.

**Note:** If necessary, place cryovials on dry ice for up to 60 minutes before thawing.

3. Immerse the cryovial in a 37°C water bath for 3 minutes (avoid submerging the cap), holding the tube stationary (no swirling). Use of a floating microcentrifuge tube rack is recommended.



Precise timing is critical to maximizing viable cell recovery.

4. Immediately remove the cryovial from the water bath, spray with 70% ethanol, and place into the biological safety cabinet.
5. Gently transfer the iCell Cardiomyocytes<sup>2</sup> cryovial contents to a sterile 50 ml centrifuge tube using a 1 ml pipettor.

**Note:** Use of a 50 ml centrifuge tube facilitates suitable mixing to minimize osmotic shock and increase cardiomyocyte viability.



Avoid repeated pipetting of the thawed iCell Cardiomyocytes<sup>2</sup> cell suspension.

6. Rinse the empty iCell Cardiomyocytes<sup>2</sup> cryovial with 1 ml of room temperature Plating Medium to recover any residual cells from the vial. Transfer the 1 ml of Plating Medium rinse from the cryovial drop-wise over 90 seconds (i.e. 1 drop every 4 - 5 seconds) to the 50 ml centrifuge tube containing the iCell Cardiomyocytes<sup>2</sup> cell suspension. Gently swirl the tube while adding the medium to mix the solution completely and minimize the osmotic shock on the thawed cells.



Drop-wise addition of Plating Medium to the cell suspension is critical to minimize osmotic shock and ensure maximum viability and attachment of the cells to the plating substrate. See the Handling iCell Cardiomyocytes Training Video available online at [www.cellulardynamics.com/cm\\_handling/](http://www.cellulardynamics.com/cm_handling/) as a reference.

## Notes

7. Slowly add 8 ml (3.5 ml for the 0.5 unit size) of room temperature Plating Medium to the 50 ml centrifuge tube. Add the first 1 ml drop-wise over 30 - 60 seconds. Then add the remaining volume over the next ~30 seconds. Gently swirl the centrifuge tube while adding the medium.



It is critical to add the Plating Medium slowly to ensure maximum viability and attachment of the cells once plated. See the Handling iCell Cardiomyocytes Training Video available online at [www.cellulardynamics.com/cm\\_handling/](http://www.cellulardynamics.com/cm_handling/) as a reference.

8. Gently mix the contents of the 50 ml centrifuge tube by inverting 2 - 3 times. Gentle mixing is critical to ensure maximum viability. Avoid vigorous shaking or vortexing of the cell suspension.

**Note:** Thaw up to 3 vials of iCell Cardiomyocytes<sup>2</sup> at one time. Once thawed, you can pool the contents of the vials before adding the rinse and final volume of Plating Medium. Follow the timing outlined in steps 6 and 7. For example, if pooling 3 vials, add each 1 ml of rinse over 90 seconds (270 seconds total).

## Chapter 6. Plating iCell Cardiomyocytes<sup>2</sup>

The recommended seeding density for iCell Cardiomyocytes<sup>2</sup> in standard cell culture plates is 156,000 viable cells/cm<sup>2</sup>. For application-specific plating instructions, see the Application Protocols available online at [www.cellulardynamics.com/lit/](http://www.cellulardynamics.com/lit/).

1. Obtain the number of viable cells/vial and viability from the Certificate of Testing.
2. Invert the thawed iCell Cardiomyocytes<sup>2</sup> cell suspension 2 - 3 times to ensure an even cardiomyocyte distribution before performing the cell count.
3. Remove a sample of cells to confirm viability using a hemocytometer (using trypan blue exclusion to identify viable cells) or an automated cell counter.
4. Calculate the final volume of room temperature Plating Medium needed to obtain a desired cell plating density (i.e.  $5 \times 10^5$  cells/ml) using the number of viable cells/vial from the Certificate of Testing. See Table 2 below for examples.

**Note:** If your application requires higher cell densities, centrifuge iCell Cardiomyocytes<sup>2</sup> at 180 x g for 5 minutes, remove the necessary amount of the Plating Medium to achieve the desired density, and gently resuspend the iCell Cardiomyocytes<sup>2</sup> pellet. Note that over-pipetting could reduce cell viability.

Culture Vessel	Surface Area (cm <sup>2</sup> )	Plating Volume (ml)	Cell Number ( $1.56 \times 10^5$ cells/cm <sup>2</sup> )
12-well Cell Culture Plate	3.8	1.2	$593 \times 10^3$
96-well Cell Culture Plate	0.32	0.1	$50 \times 10^3$

**Table 2: Summary of Recommended Volumes and Measures**

*This table provides a guide for syncytial formation only. All volumes and measures are **per well** unless otherwise indicated.*

5. Aspirate the gelatin solution from the pre-coated cell culture vessel(s).
6. Invert the cell suspension 2 - 3 times and immediately dispense the cell suspension into the pre-coated cell culture vessel(s).
7. Culture iCell Cardiomyocytes<sup>2</sup> in a cell culture incubator at 37°C, 5% CO<sub>2</sub> **for 4 hours**. Proceed immediately to Chapter 7, Maintaining iCell Cardiomyocytes<sup>2</sup>.

## Chapter 7. Maintaining iCell Cardiomyocytes<sup>2</sup>

1. Immediately before use, equilibrate iCell Cardiomyocytes Maintenance Medium (Maintenance Medium) in a 37°C water bath.
2. 4 hours post-plating iCell Cardiomyocytes<sup>2</sup>, aspirate the Plating Medium using a pipettor and replace with the appropriate volume of Maintenance Medium. Be careful not to touch or disrupt the adhered cardiomyocytes.

**Note:** Avoid excessive washing at this stage as it could affect cell adherence.

3. Replace the Maintenance Medium every other day.
4. Culture iCell Cardiomyocytes<sup>2</sup> in a cell culture incubator at 37°C, 5% CO<sub>2</sub>.

---

## Appendix A. Intellectual Property Rights, Use Restrictions, and Limited License

**A. OWNERSHIP.** The Products are covered by pending patents and patents: [www.cellulardynamics.com/patents](http://www.cellulardynamics.com/patents). Customer has a limited license to use the Products for internal research purposes for the sole benefit of the Customer, subject to the use restrictions included in subsection B of this Appendix A. Customer acknowledges and agrees that the receipt or purchase of the Products by Customer shall not be construed as a transfer of any title or the grant of any rights in or to the intellectual property embodied in the Products owned or licensed by Cellular Dynamics. In particular, no right or license to make, have made, offer to sell, or sell the Products, to modify or reproduce the Product or any part thereof, or to use the Products in combination with any other product(s), except product(s) provided or expressly licensed to Customer by Cellular Dynamics for such use, is implied or conveyed by the sale or transfer of Products to Customer.

**B. USE RESTRICTIONS.** The Products are licensed for internal research purposes only, and may not be used for any other purpose. The Products must be used in accordance with this User's Guide to which Customer, by ordering and accepting the Products, agrees. Customer shall not make, have made, offer to sell, or sell the Products. Customer shall not use the Products (or any modifications Customer makes to the Products or any cells derived, developed or expanded from the Products) in (i) the manufacture of any products, or (ii) any services for a third party. Customer may not transfer the Products (or any modifications Customer makes to the Products or any cells derived, developed or expanded from the Products) to any third party without Cellular Dynamics' prior written consent. Customer shall not reverse engineer the Products. Customer shall not use the Products, components or modifications thereof, or any cells derived, developed or expanded therefrom, in humans, in clinical trials, for diagnostic purposes involving human subjects, or for any investigational or other therapeutic use. Customer shall not use the Products directly or indirectly to derive or make any human gamete or gamete precursor cell. Customer shall use the Products in accordance with all applicable laws and regulations and any applicable institutional review board approved protocol and/or privacy office approval. Customer is not entitled to receive any data or information from Cellular Dynamics that directly identifies the donor of the biological materials from which the Products indirectly are derived or were made. Customer shall not attempt in any way to determine the identity of the donor of the biological materials from which the Products indirectly are derived or was made.

**C. DATA.** Customer agrees that if described on Customer's product quotation from Cellular Dynamics it will provide data and information as described therein to Cellular Dynamics regarding Customer's use of the Products.

---

## Appendix B. Product Provided "AS IS"

A. The Products are sold or provided "AS IS."

B. TO THE FULLEST EXTENT PERMITTED BY APPLICABLE LAW, CELLULAR DYNAMICS DISCLAIMS, ALL REPRESENTATIONS, AND WARRANTIES, EXPRESS OR IMPLIED (INCLUDING BUT NOT LIMITED TO, ANY IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE OR NONINFRINGEMENT), AND ALL LIABILITY, WITH RESPECT TO THE PRODUCTS; AND BY ORDERING AND ACCEPTING THE PRODUCTS CUSTOMER WAIVES ALL RIGHTS AND REMEDIES OTHERWISE AVAILABLE WITH RESPECT TO THE PRODUCTS.

C. Customer will be solely responsible for (i) Customer's use of the Products for a purpose or in a manner other than that for which they were designed or that is permitted or in breach of the Use Restrictions above; (ii) Customer's failure to follow this User's Guide for the use, storage, and handling of the Products however such failure is caused; (iii) Customer's failure



## Notes

to comply with any of the provisions of Appendix A above; and (iv) any abuse, other misuse or neglect of the Products by Customer or any damage or loss of the Products by events or occurrences beyond a person's (e.g., Cellular Dynamics') control including without limitation, accident, fire, vandalism and natural disasters (acts of God).

D. Customer acknowledges and agrees that Cellular Dynamics may fill Customer's order with any number of units of Products. Such units may be more units than Customer ordered. Customer will not be charged extra for any adjustments made by Cellular Dynamics. Because the number of cells in a unit may vary from lot to lot, Cellular Dynamics reserves the right to fill the order with that number of units which is sufficient to fill Customer's order.

---

### Appendix C. Limited Liability

TO THE FULLEST EXTENT PERMITTED UNDER APPLICABLE LAW, CELLULAR DYNAMICS SHALL NOT HAVE ANY LIABILITY FOR INCIDENTAL, COMPENSATORY, PUNITIVE, CONSEQUENTIAL, INDIRECT, SPECIAL OR OTHER SIMILAR DAMAGES, HOWEVER CAUSED AND REGARDLESS OF FORM OF ACTION WHETHER IN CONTRACT, TORT (INCLUDING NEGLIGENCE), STRICT PRODUCT LIABILITY OR OTHERWISE, EVEN IF CELLULAR DYNAMICS HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH DAMAGES. NOTWITHSTANDING ANY OTHER TERM OR IMPLICATION TO THE CONTRARY, UNDER NO CIRCUMSTANCES SHALL CELLULAR DYNAMICS' LIABILITY TO CUSTOMER EXCEED THE AMOUNT PAID BY CUSTOMER FOR THE PRODUCTS TO CELLULAR DYNAMICS.

© 2016 Cellular Dynamics International, Inc. All rights reserved.