


HPLC Application

Analysis of Bacitracin

Category : Drugs

Sample :
Bacitracin zinc salt

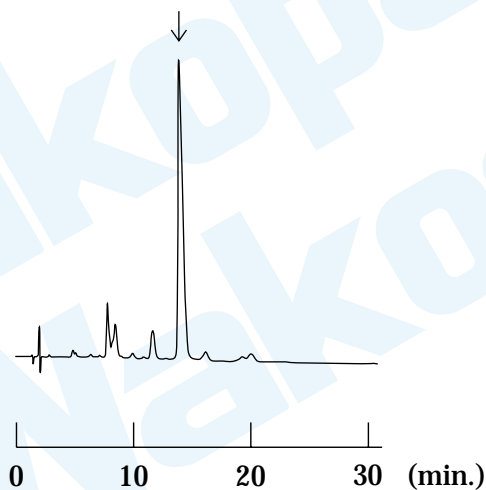

[HPLC Conditions]

Column : Wakosil-II 5C18 RS, 4.6mm I.D. x 150mm
Eluent : 20mM NaH₂PO₄, H₃PO₄ (pH 3.0)/CH₃CN = 75/25 (v/v)
Flow rate : 1.0mL/min, at 40°C
Detection : UV 210nm, 0.128AUFS.
Injection : 10 μ L, 100 μ g/mL (H₂O)